

Boule Con-5Diff A1

CONTROL

IVD


Tri-Level	72603-1K	2026/07/05
Low	72603-11	2026/07/05
Normal	72603-12	2026/07/05
High	72603-13	2026/07/05

Open vial stability: 14 days ¹

Medonic M51 / Swelab Lumi

Parameter ²	Unit ³	Low		Normal		High	
		Mean ⁴	± Range ⁵	Mean ⁴	± Range ⁵	Mean ⁴	± Range ⁵
WBC	10 ⁹ /L	3.65	± 0.50	8.27	± 1.00	18.65	± 2.50
Neu%	%	55.3	16.0	61.5	8.0	67.3	8.0
Lym%	%	36.5	21.0	29.9	8.0	20.7	12.0
Mon%	%	3.4	5.0	3.3	5.0	5.0	6.0
Eos%	%	4.8	5.0	5.3	6.0	7.0	8.0
Bas%	%	4.7	10.0	1.5	10.0	0.4	10.0
Neu#	10 ⁹ /L	2.05	0.50	5.16	0.70	12.68	2.00
Lym#	10 ⁹ /L	1.35	0.70	2.51	0.80	3.91	2.00
Mon#	10 ⁹ /L	0.13	0.30	0.28	0.50	0.94	1.10
Eos#	10 ⁹ /L	0.18	0.20	0.45	0.60	1.32	1.40
Bas#	10 ⁹ /L	0.18	0.40	0.12	0.80	0.08	1.80
RBC	10 ¹² /L	2.31	0.24	4.76	0.30	5.58	0.50
HGB	g/L	61	6	139	6	177	8
HCT	%	19.2	2.0	42.6	3.0	53.8	4.0
MCV	fL	83.2	5.0	89.6	5.0	96.5	6.0
MCH	pg	26.4	2.5	29.2	2.5	31.7	2.5
MCHC	g/L	317	30	326	30	329	30
RDW-CV	%	15.4	3.0	15.6	3.0	14.9	3.0
PLT	10 ⁹ /L	62	20	246	40	454	60
MPV	fL	8.0	3.0	7.7	3.0	8.3	3.0

No mean values and ranges are provided for research parameters-RDW-SD, PDW, PCT, P-LCR, and P-LCC. User may establish these based on their usage.

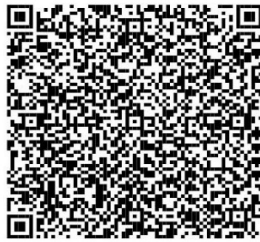

<p>1 bg Трайност на отворен флакон: 14 дни cs Stabilita otevřené ampulky: 14 dny da Åben rør stabilitet: 14 dage de Stabilität geöffneter Flaschen: 14 tage el Διατήρηση δείγματος μετά την αποσφράγιση: 14 ημέρα es Estabilidad de la cápsula abierta: 14 días et Stabiilsus avatud vialis: 14 päevade arv fr Stabilité en flacon ouvert: 14 jour hr Stabilnost otvorene bočice: 14 dani it Stabilità della fiala aperta: 14 giorni lt Stabilumas atidarius buteliukoje: 14 dieny lv Atvērt flakonū stabilitāte: 14 dienas no Åpen rør/glass stabilitet: 14 dager pl Trwałość po otwarciu fiołki: 14 Liczba dni pt Estabilidade após abertura do frasco: 14 dias ro Termenul de valabilitate al fiolei desfăcută: 14 zile ru Стабильность вскрытого флакона: 14 дней sk Stabilita otvorenej injekčnej liekovky: 14 dní sl Stabilnost odprte stekleničke: 14 dni sr Proizvod stabilan nakon otvaranja pakovanja: 14 dana sv Hållbarhet för öppnad flaska: 14 dagar tr Açık şişe stabilitesi gün: 14 gün</p>	<p>2 Параметър Parameter Parameter Parameter Παράμετρος Parámetro Parameeter Paramètre Parametar Parametro Parametras Parameters Parameter Kryterium Parâmetro Parametru Параметр Parameter Parameter Parameter Parameter Parameter</p>	<p>3 Единица Jednotka Enhed Einheit Μονάδα Unidad Ühik Unité Jedinica Unità Vienetas Vienība Enhet Jednostka Unidade Unitate Единица измерения Jednotka Enota Jedinica Enhet Birim</p>	<p>4 Средно Průměr Middelværdi Mittelwert Μέση Media Keskmine Moyenne Sredstvo Media Vidurkis Nozīmē Middelverdi Średnia Média Medie Среднее Priemer Povprečno Srednja vrednost Medelvärde Ortalama</p>	<p>5 Обхват Rozsah Område Bereich Μέση Intervalo Vahemik Intervalle Raspon Intervallo Intervalas Diapazons Område Zakres Intervalo Interval диапазон Rozsah Območje Opseg Intervall Aralığı</p>
--	--	---	--	--



www.boule.com

Boule
Boule Medical AB
Fagerstagatan 7
SE-163 53 Spånga,
Sweden

31340-13